

# The *secure* solution for Trendelenburg Positioning

Simple. Secure. Reusable.



■ **Cradles patient to maximize pressure redistribution**

■ **Unobstructed chest and IV access**

## SecureFit<sup>®</sup> TPS

The secure, easy-to-use solution that dramatically reduces setup time, resists slipping, and evenly redistributes pressure throughout the patient's body.

### Engineered to minimize slipping

The SecureFit TPS cradles the body's natural curvature, reducing peak focal pressure over the shoulder. No ancillary shoulder braces are required.

### Maximum pressure redistribution

The **TPS Gel Pad** conforms to each patient's individual anatomy and redistributes pressure to provide support while minimizing slipping.

### Simple, reusable, and easy to clean

With only a few key components and smooth, easy-to-clean surfaces, AliMed SecureFit TPS reduces setup time, which helps improve efficiency and lower costs.

### IV access simplified

**Access IV lines** quickly and easily through convenient openings between arm support and gel pad.

To order, see reverse side.

## TPS Arm Supports and Toboggans


Integrated, easy-to-clean vinyl-covered foam pads work with toboggans to protect ulnar nerve while providing easy IV access to the arms.



Toboggan

## TPS Gel Pad

Viscoelastic gel conforms to patient's anatomy and redistributes pressure for maximum support with minimal slipping.

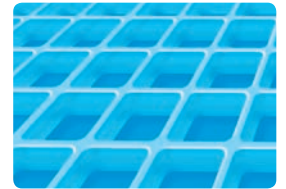

- **Simple, secure, reusable**
- **Unobstructed chest and IV access**
- **Set up in less than a minute**
- **No shoulder braces required**

## Specifications:

**TPS Foundation:** 20"W x 17"L x 6"H • **TPS Gel Pad:** 22½"W x 39"L x 1½"thick (1" gel surface with ½" Polyfoam base)  
 • Perineal cutout: 9"W x 13"L • **TPS Stabilizing Wedge Pads:** 1" Memory Foam over Polyfoam • **Cover:** Black anti-shear elastic • Antimicrobial • Fluid-proof fabric • Fluid-resistant sewn seams • CAL 117-compliant • **Stabilizing Posts:** 9¼"L x 5/8" diameter • Stainless Steel • **Toboggan:** 10"L x 8"H x 10¼"D • **Toboggan Arm Supports:** 25"L x 4½"H Polyfoam pads • Patent Pending

## Ordering Information:

#93 8040 SecureFit TPS

### Replacement Items:

- #93 8042 TPS Gel Pad
- #93 8043 TPS Stabilizing Wedge Pads, pair
- #93 8041 TPS Foundation
- #93 8038 Shoulder Wedges, pair
- #93 8044 Stabilizing Posts, pair
- #93 8045 Toboggans, pair

## TPS Stabilizing Wedge Pads

Memory foam pads minimize the incidence of poor blood flow to the shoulder.

## TPS Foundation

Strong, durable, nonskid foundation is contoured to cradle patient's natural curvature, and locks onto the surgical table with clamps for a secure hold.

