

## Worry-Free™ Alarm #71 3307

### Instructions for Application, Use, and Care

Your new AliMed® Worry-Free™ Alarm features a sealed, tamper-free battery compartment, and includes a two-year battery life guarantee. It also includes a belt clip and an adjustable strap that can be used to attach the alarm to a chair.

#### Initial Use:

When unpacking the Worry-Free™ Alarm for the first time, you will notice a small plastic tab protruding from the side of the alarm. This tab prevents the battery from draining prior to purchase. Before using the alarm, please pull this tab, then remove and discard. This activates the battery; the alarm is now on and ready to use for its 2-year life.

#### Operating Instructions:

- Connect Worry-Free™ Alarm to a bed/chair sensor pad, seatbelt, roll belt, or sensor cushion by simply plugging the connector into the jack near the bottom of the alarm.
- The system is now ready to use.
- Choose a High or Low volume tone by placing the switch in the corresponding HI or LO position.
- **Test Alarm:**
- For sensor pads: place firm pressure on the pad for at least 3 seconds. The alarm will emit two short beeps. When pressure is removed, the alarm will sound.
- For seatbelts and other similar devices: If seatbelt was already closed, open it. The alarm will sound. If seatbelt is open, close the seatbelt and you will hear the alarm emit two short beeps. When the seatbelt is subsequently opened, the alarm will sound.
- Note: If patient is placed on the sensor pad or buckled into a seatbelt, and attempts to remove the cord from the alarm, the alarm will sound.
- To reset the alarm, place pressure back on pad, close the seatbelt or similar device, or simply press the RESET button.
- Low-Battery red LED will light, along with a short beep to indicate a low battery. When this occurs, the alarm may be discarded. However, the alarm will continue to function properly under low battery for at least 5 days, allowing you time to replace it.

#### Contraindications:

- a) AliMed alarms alert caregivers to attend to patients or residents who have either risen or are attempting to rise. They do not prevent falls.
- b) Discontinue use immediately if unit is damaged or no longer functions.
- c) **Test alarm with sensor device before each use to verify that both alarm and sensor are operating properly to ensure patient safety.**
- d) **NOTICE: OPENING THE ALARM OR TAMPERING WITH BATTERIES WILL VOID THE WARRANTY.**

#### Care

- Wipe clean using standard cleaning disinfectants.
- Do not immerse in fluids.
- Check function before each use. Do not use equipment that has been damaged.

Manufactured for:

 **AliMed® inc.**  
297 High St., Dedham, MA 02026 USA  
PH 800-225-2610 • FAX 800-437-2966

AliMed.com  
#7140REV1213 CCO681  
©2013 AliMed® inc. All rights reserved.  
Made in China

## Worry-Free™ Alarm #71 3307

### Instructions for Application, Use, and Care

Your new AliMed® Worry-Free™ Alarm features a sealed, tamper-free battery compartment, and includes a two-year battery life guarantee. It also includes a belt clip and an adjustable strap that can be used to attach the alarm to a chair.

#### Initial Use:

When unpacking the Worry-Free™ Alarm for the first time, you will notice a small plastic tab protruding from the side of the alarm. This tab prevents the battery from draining prior to purchase. Before using the alarm, please pull this tab, then remove and discard. This activates the battery; the alarm is now on and ready to use for its 2-year life.

#### Operating Instructions:

- Connect Worry-Free™ Alarm to a bed/chair sensor pad, seatbelt, roll belt, or sensor cushion by simply plugging the connector into the jack near the bottom of the alarm.
- The system is now ready to use.
- Choose a High or Low volume tone by placing the switch in the corresponding HI or LO position.
- **Test Alarm:**
- For sensor pads: place firm pressure on the pad for at least 3 seconds. The alarm will emit two short beeps. When pressure is removed, the alarm will sound.
- For seatbelts and other similar devices: If seatbelt was already closed, open it. The alarm will sound. If seatbelt is open, close the seatbelt and you will hear the alarm emit two short beeps. When the seatbelt is subsequently opened, the alarm will sound.
- Note: If patient is placed on the sensor pad or buckled into a seatbelt, and attempts to remove the cord from the alarm, the alarm will sound.
- To reset the alarm, place pressure back on pad, close the seatbelt or similar device, or simply press the RESET button.
- Low-Battery red LED will light, along with a short beep to indicate a low battery. When this occurs, the alarm may be discarded. However, the alarm will continue to function properly under low battery for at least 5 days, allowing you time to replace it.

#### Contraindications:

- a) AliMed alarms alert caregivers to attend to patients or residents who have either risen or are attempting to rise. They do not prevent falls.
- b) Discontinue use immediately if unit is damaged or no longer functions.
- c) **Test alarm with sensor device before each use to verify that both alarm and sensor are operating properly to ensure patient safety.**
- d) **NOTICE: OPENING THE ALARM OR TAMPERING WITH BATTERIES WILL VOID THE WARRANTY.**

#### Care

- Wipe clean using standard cleaning disinfectants.
- Do not immerse in fluids.
- Check function before each use. Do not use equipment that has been damaged.

Manufactured for:

 **AliMed® inc.**  
297 High St., Dedham, MA 02026 USA  
PH 800-225-2610 • FAX 800-437-2966

AliMed.com  
#7140REV1213 CCO681  
©2013 AliMed® inc. All rights reserved.  
Made in China

## Worry-Free™ Alarm #71 3307

### Instructions for Application, Use, and Care

#### To Reorder:


- #71 3307 Worry-Free™ Alarm unit only
- #71 330710 Worry-Free™ Alarm, 10 pack
- #71 3310 Worry-Free™ Alarm with 2 year chair sensor pad
- #71 3311 Worry-Free™ Alarm with 2 year bed sensor pad
- #71 3312 Worry-Free™ Alarm with E-Z release seatbelt
- #71 3313 Worry-Free™ Alarm with Early Warning E-Z release seatbelt
- #71 3314 Worry-Free™ Alarm with Buckled Seatbelt
- #71 3315 Worry-Free™ Alarm with Breakaway Lap Cushion
- #71 3316 Worry-Free™ Alarm with Roll Control Belt
- #2032 VOICE Alarm™

### Pull-cord and sensor alarm in one!

**VOICE Alarm™** allows caregivers to record a personalized message which plays when the alarm is triggered. Perfect choice for patients with dementia, those who are easily confused or disoriented by loud sirens, or patients who need instructions in a foreign language. Uses a magnetic pull-cord or connects to a sensor pad—your choice! Select by flipping a switch. You also get the option of selecting a regular alarm tone, with two volumes. Includes mounting holster.

85- or 100-dB alarm setting • 20-second recordable alarm alert • Four different alerts • Silent feature when used in combination with the nurse-call jack alarm; can be set to silent in patient's room • Low-battery light • Mounting holster • Nurse-call jack • For use as a pull-cord or with sensor devices—sensor pads, bed/chair, seatbelts, bed belt, cushions • Adjustable pull-cord • 4 AA batteries, included

#2032 VOICE Alarm™

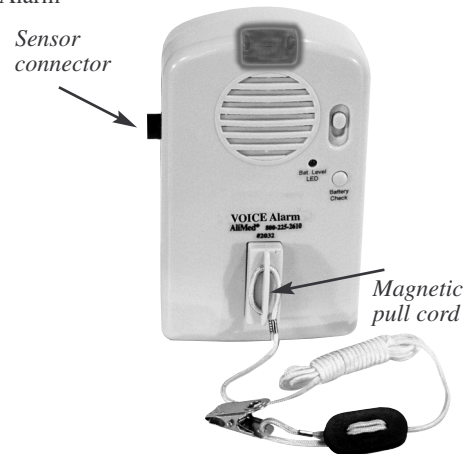

Manufactured for:

 **AliMed® inc.**  
297 High St., Dedham, MA 02026 USA  
PH 800-225-2610 • FAX 800-437-2966

AliMed.com  
#7140REV1213 CCO681  
©2013 AliMed® inc. All rights reserved.  
Made in China

## Worry-Free™ Alarm #71 3307

### Instructions for Application, Use, and Care

#### To Reorder:

- #71 3307 Worry-Free™ Alarm unit only
- #71 330710 Worry-Free™ Alarm, 10 pack
- #71 3310 Worry-Free™ Alarm with 2 year chair sensor pad
- #71 3311 Worry-Free™ Alarm with 2 year bed sensor pad
- #71 3312 Worry-Free™ Alarm with E-Z release seatbelt
- #71 3313 Worry-Free™ Alarm with Early Warning E-Z release seatbelt
- #71 3314 Worry-Free™ Alarm with Buckled Seatbelt
- #71 3315 Worry-Free™ Alarm with Breakaway Lap Cushion
- #71 3316 Worry-Free™ Alarm with Roll Control Belt
- #2032 VOICE Alarm™

### Pull-cord and sensor alarm in one!

**VOICE Alarm™** allows caregivers to record a personalized message which plays when the alarm is triggered. Perfect choice for patients with dementia, those who are easily confused or disoriented by loud sirens, or patients who need instructions in a foreign language. Uses a magnetic pull-cord or connects to a sensor pad—your choice! Select by flipping a switch. You also get the option of selecting a regular alarm tone, with two volumes. Includes mounting holster.

85- or 100-dB alarm setting • 20-second recordable alarm alert • Four different alerts • Silent feature when used in combination with the nurse-call jack alarm; can be set to silent in patient's room • Low-battery light • Mounting holster • Nurse-call jack • For use as a pull-cord or with sensor devices—sensor pads, bed/chair, seatbelts, bed belt, cushions • Adjustable pull-cord • 4 AA batteries, included

#2032 VOICE Alarm™



Manufactured for:

 **AliMed® inc.**  
297 High St., Dedham, MA 02026 USA  
PH 800-225-2610 • FAX 800-437-2966

AliMed.com  
#7140REV1213 CCO681  
©2013 AliMed® inc. All rights reserved.  
Made in China